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SPA PRO

Special Jade massage beds

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Special model s MP3 prehrávačom

Special modell with MP3 player

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Table of contents:

Product Introduction	
Main Features	04
Messages For Users	04
Contraindications	04
Precautions	05
Inspection Before Using The Product	05
Inspection After Using The Product	05
Remote Controller	06
Instructions For Using LCD Controller	07
Program Diagram	08, 09, 10, 11
How To Use Our Product	
Using The Product	12
Installation Of The Product In Your Home	
Points To Remember Before Installation	12
The Set-up Procedure	13
Replace the fuse	13
Maintenance & Cautions	
Safety Practice	14
Care And Storage	14
Remote Controller Operation & Error Messages Explanations	
Service	15
Specifications	16

Product Introduction

Main features

1. Four or nine groups of premium jade internal projector for maximum effect
2. Two external jade projectors with 5/9 balls to be applied on the parts of body unreachable by the bed
3. Back massage mattress can be adjusted to 45 degree angle to suit for your personal needs
4. Heated jade roller reach temperatures from 40 to 70 to effectively transfer heat to your needed discomfort spots on your back
5. Neck jade roller with linear rail to provide extra treatment for the head and neck discomfort
6. Tourmaline stone has five major benefit for human body
 - a. Release of 4-14um high energy far infrared ray
 - b. Release of negative ion permanently
 - c. Generate 0.06ma bio-current
 - d. Rich in minerals and elements
 - e. Eliminates body odor and kills bad bacteria.
7. Chiropractic, acupressure and infrared therapy
8. With MP3 music function.

Messages for users

Used as supplementary instrument treatment for lumbar spondylosis (chronically strain of lumbar muscles, slipped lumbar disc, back pain) and cervical spine problems

Contraindications

People with the follow illness are strictly prohibited to use this product:

1. Pregnant women
2. With tumor
3. Severe osteoporosis
4. Severe high blood pressure
5. Acute suppurative inflammation
6. Hemorrhagic tendency
7. Artificial cardiac pacemaker
8. Bone fracture
9. Hypersensitive to electricity heat

Precautions

1. Cotton clothes are recommended during therapy. Please remove the belt and tie before therapy
2. Please do not place external projectors directly on the skin to avoid the possible burn in case of excessive heat resulted by improper setting of the temperature.
3. Please do not place external projectors directly over your face to avoid the contact with eyes.
4. If you are sensitive to heat, please lower the temperature to proper level before using the product.
5. Please follow the instructions and do not use other brand products while using this products.
6. Please do not go to sleep during the program.
7. This device can not replace any other medical treatments/surgeries suggested or prescribed by your doctor for curing/healing purpose. If you symptom persists or worsen, please seek for medical help.


Inspection before using the product

1. Make sure your power cable is plugged in a working power outlet;
2. Check if external projectors are plugged in properly with the bed;
3. Check if main power switch and remote controller work properly;
4. Check if temperature display is correct on the remote controller;
5. Check if all the lamps are working properly;
6. Please make sure the projectors temperature is around 45 -55 ;
7. Check if any abnormality about projectors;

Inspection after using the product

1. Check if temperature display is turned off;
2. Check if all lamps are shut off;
3. Check if power switch is at off position;
4. Check if power cable is unplugged from power outlet;

Remote Controller

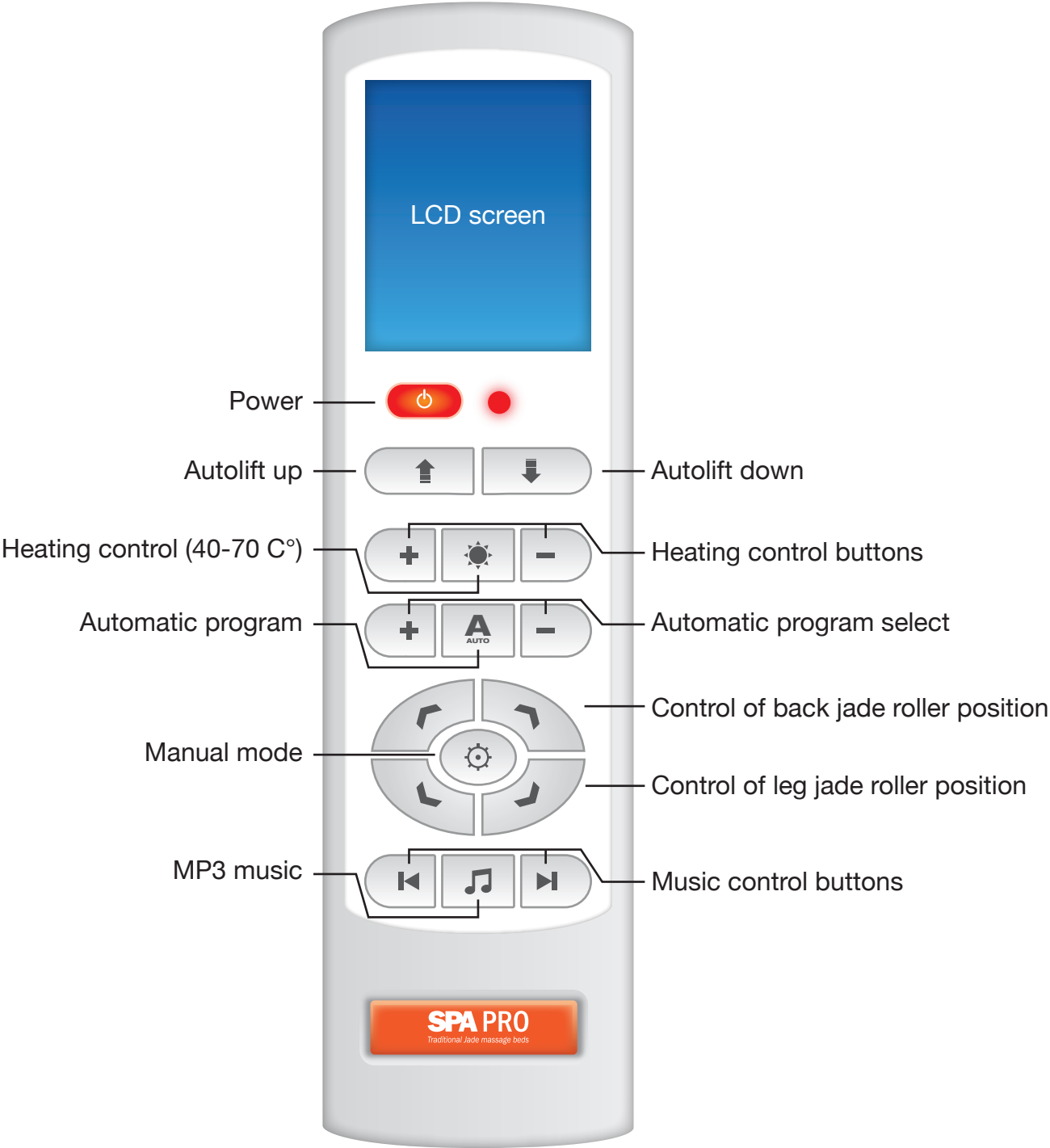Instructions for using LCD controller

Power button:
press this button, the bed enters working state, the default working model is model 1 and default working time is 40mins.

Autolift up button:
press this button, the bed back mattress would be up if bed is working state.

Autolift down button:
press this button, the bed back mattress would be down if bed is working state.


Temperature button:
Press this button continuously, can select the different heating device, no operation within 5 seconds, the heating device is confirmed, the default temperature is 55 , press the left button” +” or the right button” -“,increasing or decreasing the temperature from 40 to 70 , press one time, the temperature changes 5 .

Auto model button:
Press this button, can select the different working model N1,N2,N3,N4,N5,N6,N7, the default model is N1, press the left button” +” or the right button” -“,select different model, no operation within 5 seconds, the model is confirmed, then press”+” or”-“ button, can change the working time,press one time, the time changes one minute, the default time is 40mins, the longest time is 99mins.

Manual model button:
Press this button, the bed enters manual model state.

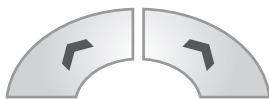Moving jade roller button:
Under manual model state, press this button, can move the back or leg jade roller up.

Moving jade roller button:
Under manual model state, press this button, can move the back or leg jade roller down.

Mp3 music button:
Press this button, the music start working, press again, the Mp3 stop working.


Music selecting button:
Press this button, select the front music.


Music selecting button:
Press this button, select the back music.

Program:

1 Program:

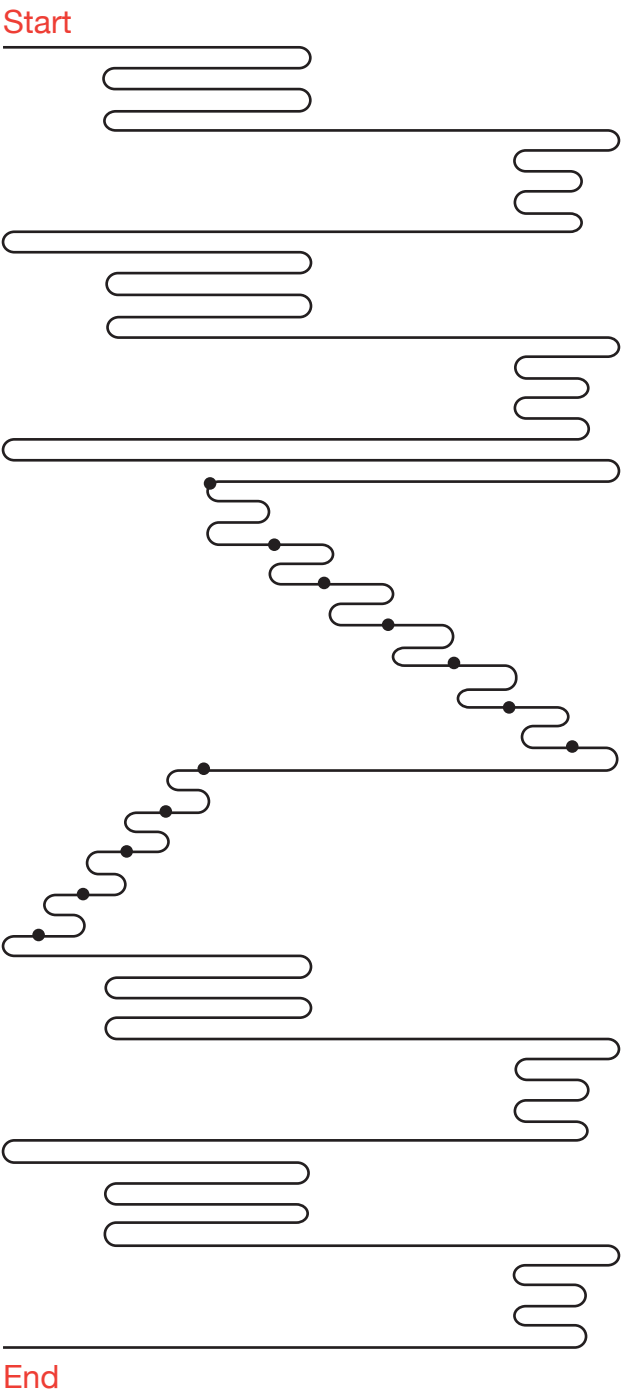
6. Program:

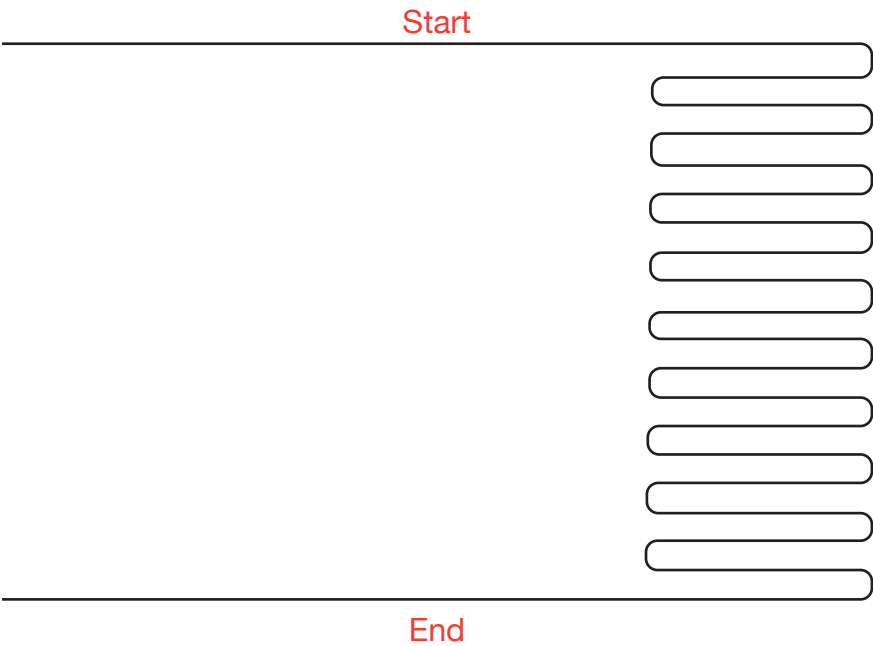How to use our product


- a. See the picture below. Lay with your belly up on the massage bed. Keep your head fist down the top of the massage bed. Rest your legs and arms on the mattress.
- b. When using Automatic modes, please place the external projectors on your desiredn spots (Refer to below description for the safety practice)
- c. You may continuously use Automatic or manual mode to reach the desired results.

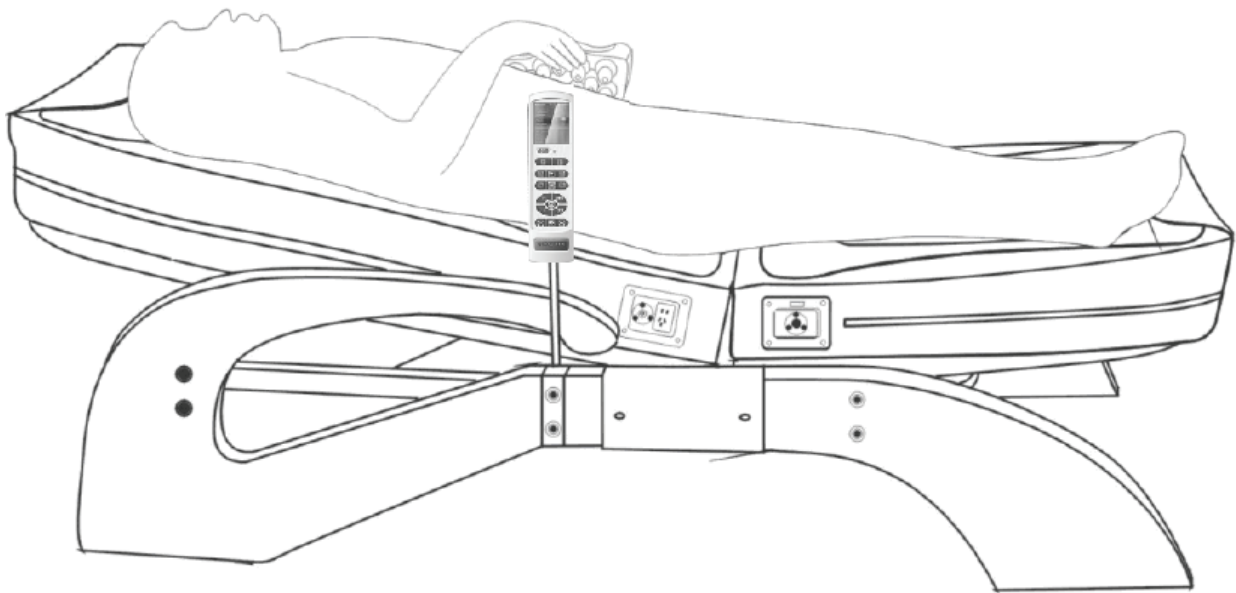Installation Of The Product In Your Home

- a. See the picture below. Lay with your belly up on the massage bed. Keep your head fist down the top of the massage bed. Rest your legs and arms on the mattress.
- b. When using Automatic modes, please place the external projectors on your desiredn spots (Refer to below description for the safety practice)
- c. You may continuously use Automatic or manual mode to reach the desired results.

The Set-up Procedure

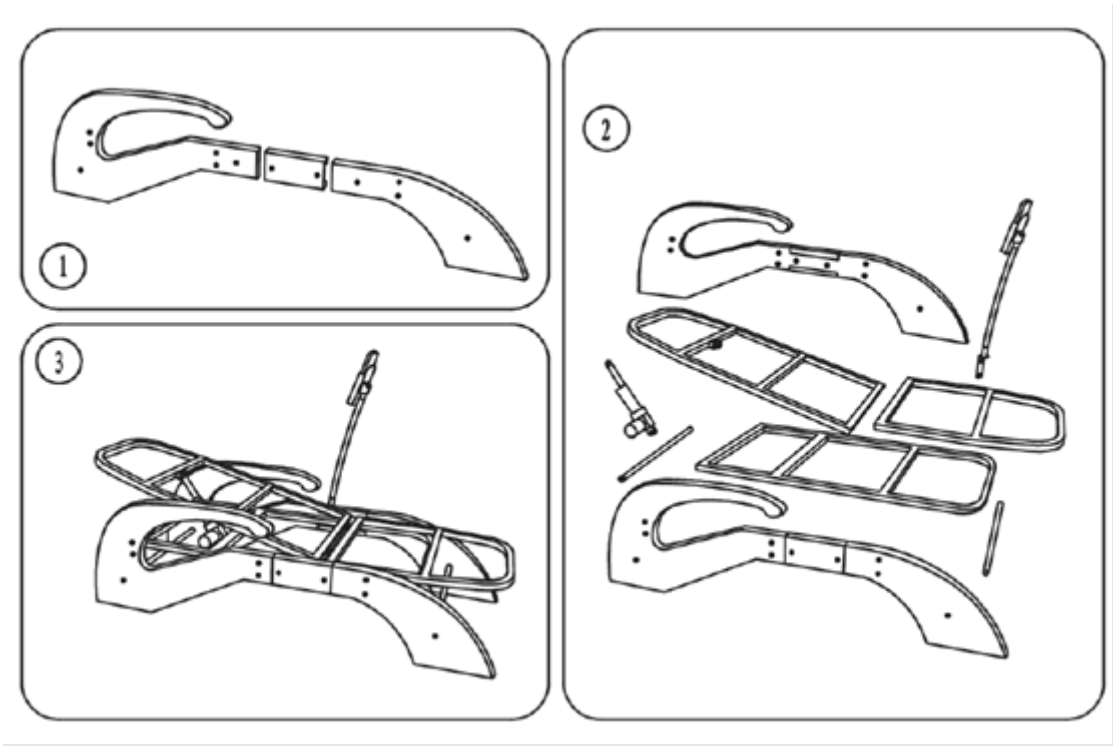
Replace the fuse

1. Draw out the power cable
2. Ascertain the upper covered trough in the recess,then pull out the socket of fuse by screw flat tip screw driver

3. Take out the spare fuse gently and replace it.

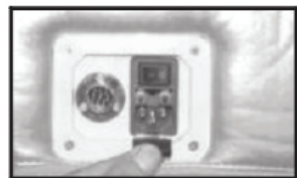4. Insert the fuse in the clutch base according to the primary direction.

Maintenance & Cautions

Safety practice

1. This product has been designed for your maximum comfort and safety. However, improper use may result in fire and electricity shock.
2. Please make sure your power output is AC 120 V/220V before each usage.
3. Please use power cable that shipped with the product. Do not use other power cables other than supplied.
4. Please go to your local store or your service provider to replace the power cable if it is damaged.
5. Do not overload your power.
6. Please do not urzip or open the leather cover in case of danger when the bed is in idle or operation.

Care and storage

1. Please unplug the power cable after each usage.
2. Do not attempt to open the product if you are not a professional service engineer
3. In the case of abnormal noises, please stop using the production immediately and contact the service people.
4. Please do not mix use all cables to avoid danger.
5. Please do not pull directly on power cable when it is difficult to unplug to avoid danger.
6. Do not use wet hands to touch cables to avoid electricity shock.
7. Use separate outlet with other electronics.
8. When idle for a long period of time or in the storm weather, please unplug the power cables from your power outlet.

Remote controller operation & error messages explanations

Safety practice

The concrete question please examine the maintenance instructions handbook.

Self troubleshooting

Problems	Check	Solution
Remote control has no or wrong display	Power cable plug in properly	Reconnect power cable to out let
	Remote control plug in properly	Reconnect remote control cable to the outlet
	Burnt fuse	Replace the fuse
	Power switch on	Turn power switch to On position
	Press wrong Function buttons	Read the instruction
Projector does not light up	Check if power is plugged in and outlet has power	Call service if none of these work. Do not open the product yourself
Internal projector does not move	Press the remote control to choose the programs	Shut off power and restart the product
	No indicator light	Shut off power and restart the product

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